


RANGE 3/5 NF EN ISO 1182



SAFETY STORAGE CABINETS

for Toxic, Noxious or Flammable products

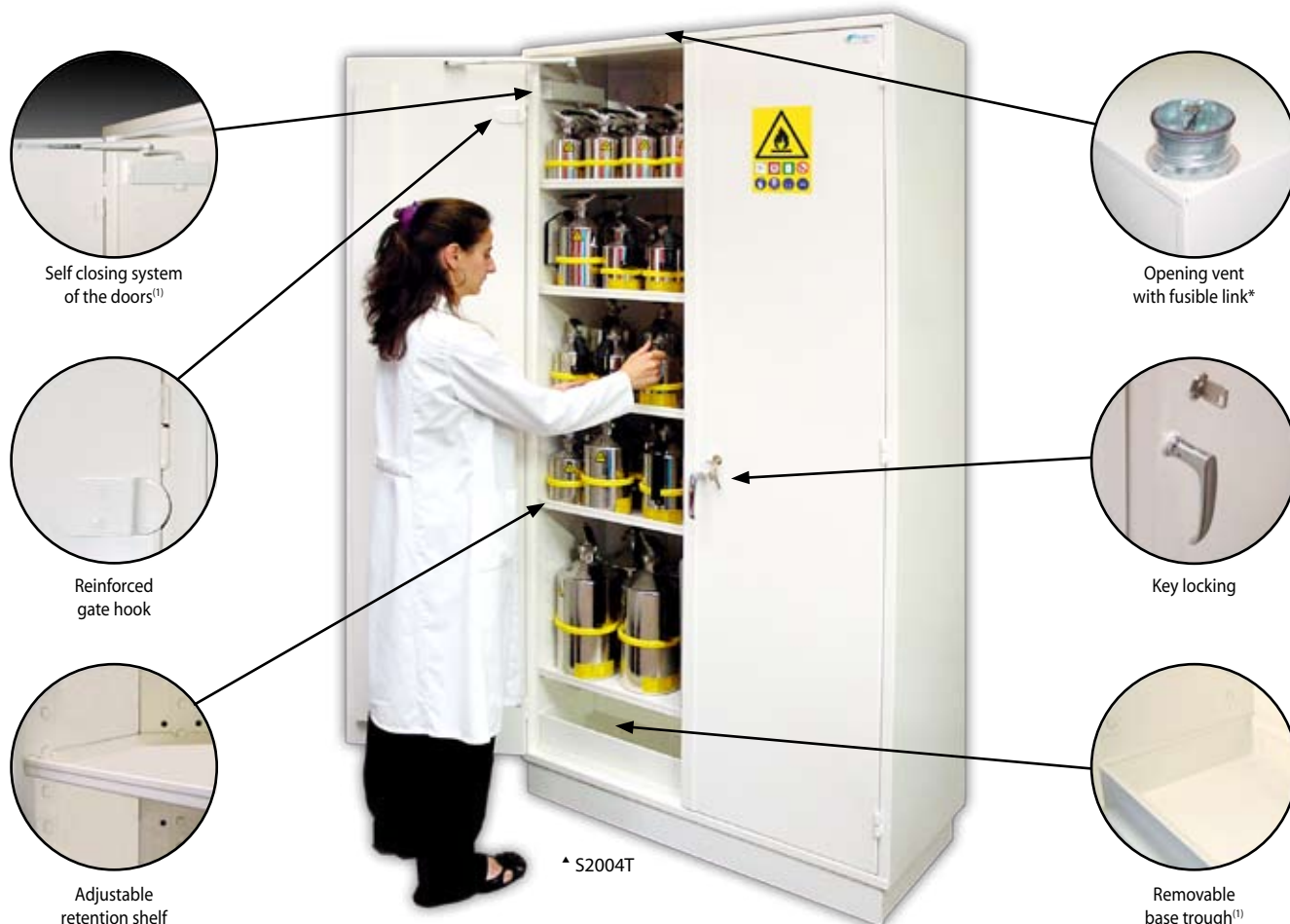

Self closing system of the doors⁽¹⁾

Opening vent with fusible link*

Reinforced gate hook

Key locking

Adjustable retention shelf

Removable base trough⁽¹⁾

^ S2004T

DIRECTIVES & NORMS

- In as per EN ISO 1182.
- Normalized labelling according to ISO 3864 norm and 92/58/CEE directive.

ADVANTAGES

- Large storage capacity.
- Cabinets for various families of products storage (noxious, toxics and flammables) with same design for a better homogeneity.

PASSIVE SECURITY

- Built in double wall steel 12/10th.
- 20 mm thickness of mineral fibre insulation.
- White epoxy paint RAL 9010.
- «Toxics», «Noxious» or «Flammables» large warning label.
- Adjustable containment shelves.

ACTIVE SECURITY

- Inside lateral ventilation shaft with high opening and connecting piece (Ø 100 mm).
- Ventilation outlet on the roof with connecting sleeve (Ø 100 mm).
- Reinforced gate hook, sealing off by 2 points lock.
- Key locking.

OPTIONS

- Additional retention shelf (G03, G04 or G06).
- Containment shelf in PVC (EP03, EP04 or EP06).
- Base through in steel (GL03, GL04, or GL06).
- Base through in PVC (BP03, BP04 or BP06).
- Absorbing carpets (ref. JF10).



ADDITIONAL EQUIPMENT OF THE RANGE 5

- Removable base through in steel at the bottom.
- Self closing system of the doors.
- Opening vent with fusible link.

RANGE 5 OPTIONS

- Open doors locking system with self closing system at 70°C.
- Additional retention shelf (G03, G04 ou G06).
- Containment shelf in PVC (EP03, EP04 ou EP06).
- Containment sump in PVC (BP03, BP04 ou BP06).


Ventilation and Filtration Boxes

To protect the users health and to improve the life cycle of your safety cabinet, we recommend to install an exhausting system. Do not hesitate to contact-us for more information.



RANGE 3/5

NF EN ISO 1182

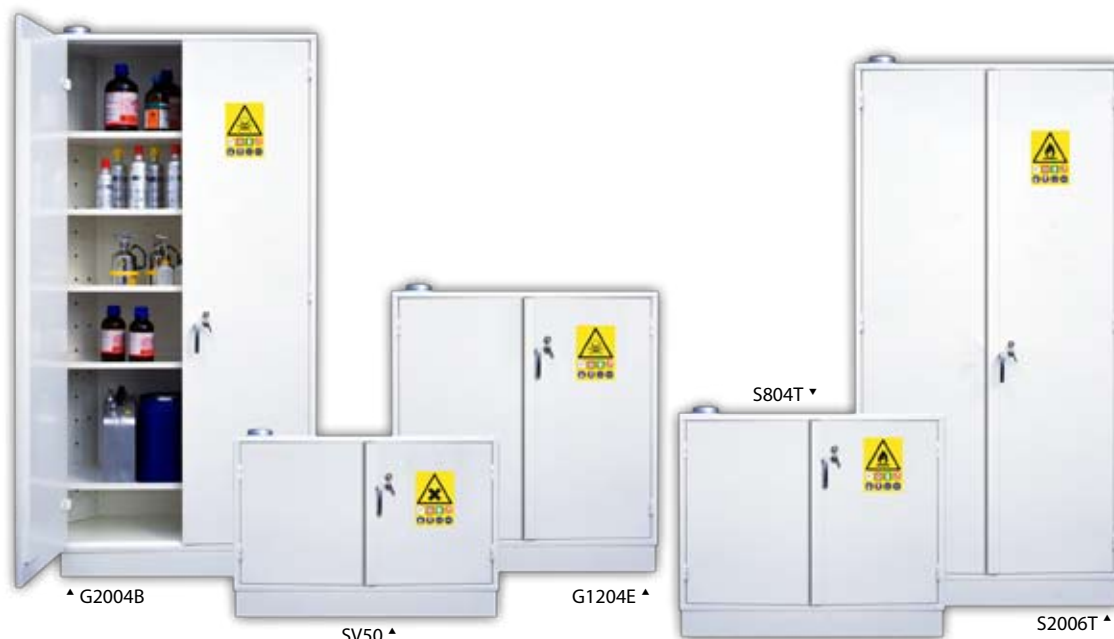

Flammable cabinets

Corrosive cabinets

Toxic cabinets and multirisk

Pesticide range

Fume hoods, filtering cabinets - Ventilation

Containment and Cans

File cabinets Anti-fire equipments

Showers and first-aid equipments

The manufacturer reserves the right to discontinue or change specifications or designs at any time without notice and without incurring obligations. 05/2017

RANGE 3

P/N	Description	Outside dim. H x W x D (mm)	Inside dim. H x W x D (mm)	Storage capacity (L)	Number of shelves	Maximum load by shelf ⁽¹⁾ (kg)	W. (kg)
SV50	Under-bench cabinet (depth: 450 mm)	750 ⁽²⁾ x 1000 x 450	650 x 950 x 400	50	1	75	62
G803G	Low cabinet (depth: 300 mm)	895 x 1000 x 300	745 x 950 x 250	70	1	75	58
G804G	Low cabinet (depth: 450 mm)	895 x 1000 x 450	745 x 950 x 400	100	1	75	68
G806G	Low cabinet (depth: 600 mm)	895 x 1000 x 600	745 x 950 x 550	135	1	75	82
G1204E	Working cover cabinet (depth:450mm)	1295 x 1000 x 450	1145 x 950 x 400	150	3	75	92
G2004B	Tall cabinet (depth: 450 mm)	2095 x 1000 x 450	1945 x 950 x 400	250	5	75	145
G2006B	Tall cabinet (depth: 600 mm)	2095 x 1000 x 600	1945 x 950 x 550	300	5	75	166

RANGE 5

P/N	Description	Outside dim. H x W x D (mm)	Inside dim. H x W x D (mm)	Storage capacity (L)	Number of shelves	Maximum load by shelf ⁽¹⁾ (kg)	W. (kg)
S604T	Under-bench cabinet (depth: 450 mm)	750 ⁽²⁾ x 1000 x 450	650 x 950 x 400	50	1	75	69
S804T	Low cabinet (depth: 450 mm)	895 x 1000 x 450	745 x 950 x 400	100	1	75	75
S806T	Low cabinet (depth: 600 mm)	895 x 1000 x 600	745 x 950 x 550	135	1	75	90
S1204T	Working cover cabinet (depth:450mm)	1295 x 1000 x 450	1145 x 950 x 400	150	3	75	99
S2004T	Tall cabinet (depth: 450 mm)	2095 x 1000 x 450	1945 x 950 x 400	250	5	75	152
S2006T	Tall cabinet (depth: 600 mm)	2095 x 1000 x 600	1945 x 950 x 550	300	5	75	174

OPTIONS

P/N	Description	W. (kg)
G03	Additional shelf for G803G	3
G04	Additional shelf for SV50-G804G-G1204GE-G2004B-S604T-S804T-S1204T-S2004T	4
G06	Additional shelf for G806G-G2006B-S806T-S2006T	5
EP03	Containment shelf in PVC for G803G	6
EP04	Containment shelf in PVC for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	6
EP06	Containment shelf in PVC for G806G-G2006B-S806T-S2006T	6
GL03	Containment sump for G803G	4
GL04	Containment sump for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	5
GL06	Containment sump for G806G-G2006B-S806T-S2006T	6
BP03	Containment sump in PVC for G803G	3
BP04	Containment sump in PVC for SV50-G804G-G1204E-G2004B-S604T-S1204T-S2004T	4
BP06	Containment sump in PVC for G806G-G2006B-S806T-S2006T	5
FTR5	Self closing system at 70°C with open doors locking system (must be ordered together with the cabinet)	-

⁽¹⁾ Evenly distributed ⁽²⁾ In option: height 670 mm. Contact-us for more information.